

Invia[®] Motion[™] NPWT System

PERSONAL-PUMP CONVENIENCE WITH REUSABLE PUMP PERFORMANCE

Personal-pump convenience

- **Compact** design for **comfort** and **mobility**
- **Simple** inventory management
- **No complex rental fees**

Reusable pump performance

- **Intelligent Pressure Control** delivers set pressure at the wound site¹
- **Dynamic Exudate Removal** helps prevent blockages¹
- **Canister based system** for efficient exudate management




Invia Motion provides an efficient solution to support clinicians with Negative Pressure Wound Therapy management in all patient care settings.


The Invia Motion delivers the convenience of a personal pump while offering clinical flexibility with a choice of pressure settings, therapy modes and compatibility with both foam and gauze dressings. The Invia Motion is appropriate for use in all patient care settings. The compact design and quiet operation helps patients to comfortably go about daily activities for the duration of therapy.

Invia Motion delivers therapy confidence. With a double lumen tubing design from the pump to the wound, the pump actively controls the set pressure at the wound site. In addition, small blockages in the tubing are resolved by the pump with assurance that clinicians and patients will see and hear notifications for any unresolved blockages.


The **Quick-connector** interface provides **easy and secure** attachment² between the dressing tubing and the canister tubing.

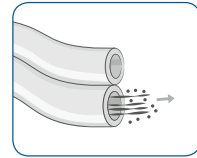

Double lumen to the wound

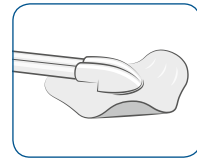

FitPad

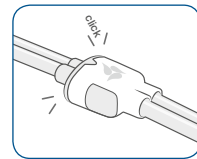

Quick-connector

Key features at a glance

- Providing NPWT in the hospital, in transition, and at home
- Clear user interface supported by 4 buttons and a digital display
- Discreet carrying case, battery autonomy of 10 hours and compact design for comfort and mobility
- Preset pressure settings from -40 to -175mmHg (-125mmHg default setting)
- Easy and secure tubing attachment with the Quick-connector
- Choice of therapy runtimes to meet clinical needs



Technical data



Medium vacuum
- 40 mmHg to -175 mmHg



Battery life
> 10 Hours

Modes: Constant and Intermittent
Tubing: Double Lumen Tubing
Canister: 150 ml Canister/Tubing set



Low flow
1 l/min.



350 g
0.88 lbs
without canister



H x W x D
99 x 98 x 52 mm
3.90 x 3.86 x 2.05 inch

2-year warranty

Find out more about Medela NPWT portfolio at www.medela-healthcare.us

CAUTION: U.S. Federal law restricts this device to sale by or on the order of a physician.

¹ Performance data on file at Medela AG
² Usability Validation data on file at Medela AG

Local contact:

Medical Vacuum Technology
for Healthcare Professionals

Please contact us or your local Medela representative for details.

Medela and Invia are registered in the US Patent and Trademark Office. Invia Motion is a trademark of Medela Holding AG, Switzerland. No trademark registrations in Canada.

 Medela AG
Lättichstrasse 4b
6341 Baar, Switzerland
www.medela.com

USA
Medela LLC
1101 Corporate Drive
McHenry, Illinois 60050
USA
Phone +1 877 735 1626
Fax +1 815 307 8942
info-healthcare@medela.com
www.medela-healthcare.us

Canada
Medela Canada Inc.
4160 Sladeview Cres., Unit #8
Mississauga, Ontario, L5L 0A1
Canada
Phone +1 800 435 8316
Fax +1 800 995 7867
info@medela.ca
www.medela.ca